

# **P515/535**

## **Emitter-Detector Module**

### **User's Instruction**

1. Edition: March 17, 2011  
515\_535\_module\_01.doc

© Heinz Walz GmbH, 2011


Heinz Walz GmbH • Eichenring 6 • 91090 Effeltrich • Germany  
Phone +49-(0)9133/7765-0 • Telefax ++49-(0)9133/5395  
E-mail [info@walz.com](mailto:info@walz.com) • Internet [www.walz.com](http://www.walz.com)



# P515/535

## Emitter-Detector Module

for simultaneous measurements of 515-520 nm ("P515")  
and 535 nm ("scattering") absorbance changes

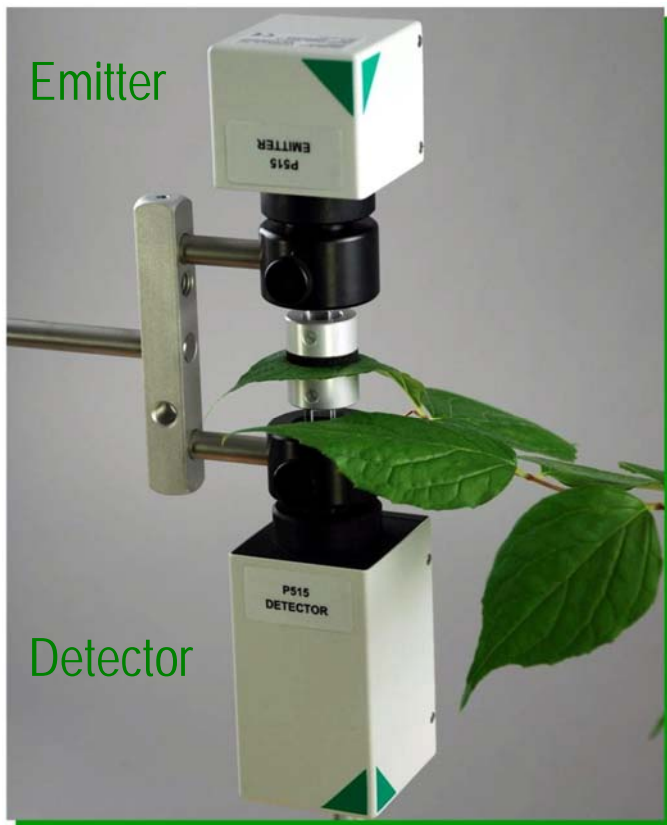



# CONTENT

1	P515/535 Emitter-Detector Module	1
1.1	Extent of Delivery	2
1.1.1	P515/535 Emitter Head (DUAL-EP515)	2
1.1.2	P515/535 Detector Head (DUAL-DP515)	2
1.2	Setup for parallel assessment of P515 and 535 nm scattering	4
1.3	Setup for parallel assessment of P515 and chlorophyll fluorescence	6
1.4	Spectra of measuring light, P515 and 535 nm, scattering changes	8
1.4.1	Data source for Figure 3	10
2	Accessory: DUAL-K25 Quartz Glass Cuvette for Suspensions cuvette	11
3	Index	13

## APPENDIX

4	Schreiber U, Klughammer C (2008) New accessory for the Dual-PAM-100: The P515/535 module and examples of its application. PAN 1: 1-10	
---	---	--



# 1 P515/535 Emitter-Detector Module

This manual provides a brief introduction to the use of P515/535 Emitter-Detector Module. More detailed information is given in the Appendix of this booklet by highly informative article written by Ulrich Schreiber and Christof Klughammer (PAN 2008, 1: 1 - 10, URL:

[http://www.walz.com/products/chl\\_p700/pan/overview.html](http://www.walz.com/products/chl_p700/pan/overview.html)).

Generally speaking, the P515/535 emitter-detector module is designed to measure absorbance changes in the green spectral region (500 - 550 nm) of photosynthetic samples.

In this range, the two most prominent absorbance changes originating from primary photosynthetic processes are the P515 electrochromic band shift and a 535 nm scattering change. The P515 band shift is induced by changes in the electric field across photosynthetic membranes, and the 535 nm scattering change arises from conformational variations during acclimation of photosynthesis to light. In addition, long-term measurements of P515 can be affected by absorbance changes caused by the xanthophyll cycle.

In the green spectral region, the maximum P515 absorbance change occurs between 515 and 520 nm; the maximum apparent absorbance difference caused by the scattering change is situated around 535 nm (compare Fig. 3B and C). In comparison, the maximum absorbance difference caused by carotenoid conversions within the xanthophyll cycle is at 505 nm.

The P515/535 emitter-detector module measures the P515 in the dual wavelength mode with 520 nm as the sample and 550 nm as the reference wavelength. The 535 nm scattering change is measured in the single wavelength mode as variation in attenuation of a 535 nm beam of measuring light. For all

the measuring light sources, typical intensity spectra are depicted in Fig. 3.

## 1.1 Extent of Delivery

### 1.1.1 P515/535 Emitter Head (DUAL-EP515)

**Green measuring light:** 520 nm sample wavelength and 550 nm reference wavelength for dual wavelength P515 measurements, 535 nm for single wavelength measurement of 535 nm scattering changes.

**Actinic light:** Chip-on-board (COB) LED array emitting maximally at 635 nm and delivering up to  $2\,500\ \mu\text{mol m}^{-2}\ \text{s}^{-1}$  of photosynthetically active radiation (PAR), single turnover flashes up to an intensity of  $200\,000\ \mu\text{mol m}^{-2}\ \text{s}^{-1}$  PAR and adjustable between 5 and 50  $\mu\text{s}$ , and multiple turnover flashes up to  $20\,000\ \mu\text{mol m}^{-2}\ \text{s}^{-1}$  PAR and adjustable between 1 and 1000 ms. Far-red LED emitting maximally at 740 nm for preferential excitation of PS I.

**Dimensions:** 11 cm x 5.0 cm x 5.5 cm (L x W x H).

**Weight:** 400 g (incl. cables, 1 m long).

### 1.1.2 P515/535 Detector Head (DUAL-DP515)

**Signal detection:** 10 x 10 mm PIN-photodiode equipped with a 2 mm blue-green filter sandwich (BG39 and DT-Cyan) transmitting all three measuring light wavelengths.

**Dimensions:** 15.5 cm x 5.0 cm x 5.5 cm (L x W x H).

**Weight:** 350 g (incl. cable, 1 m long).



### 1.1.3 Absorption Foil (“Artificial Leaf”)

2.0 cm x 2.4 cm rectangle of Roscolux #01 - Light Bastard Amber color foil (Rosco Laboratories Inc., Stamford, CT, USA) wrapped in light-scattering paper.

## 1.2 Setup for parallel assessment of P515 and 535 nm scattering

Evaluation of P515 and the 535 nm scattering changes requires that the transmission of measuring light by the sample is measured. For these light transmission measurements, the measuring light source (DUAL-EP515), the sample, and the measuring light detector (DUAL-DP515) must form a line as Fig. 1 schematically outlines (see Setup and Connections).

Also, Fig. 1 depicts how emitter and detector are connected to the DUAL-C Power-and-Control-Unit: the DUAL-EP515 emitter is plugged into the 6 pin socket for measuring light (EMITTER) and power for actinic light is provided by the lower LED array socket (LED ARRAY, socket E). The DUAL-DP515 photodiode detector is plugged into the socket DETECTOR 1.

The DualPAM software provides a file denoted Walz\_P515.DEF which contains instrument settings optimized for P515 measurements. These settings can be activated using the button “Open User Settings” followed by opening the Walz\_P515.DEF file.

For parallel measurements of P515 and 535 nm scattering, select DP515/535 as Detector 1 in the mode-setting window of the DualPAM software. Check N.C. (none connected) in the Detector 2 section (see Fig. 1). Finally, select I(535) + P515 as the Dual Channel mode.

## Mode-setting window of DualPAM software

<b>Measure Mode</b>		<b>Detector Type</b>	
Single Channel	Dual Channel	Detector 1	Detector 2
<input type="radio"/> P515	<input checked="" type="radio"/> I(535nm) + P515	<input type="radio"/> DB	<input checked="" type="radio"/> N.C.
<input type="radio"/> I(535)	<input type="radio"/> P515 + Ext 1	<input type="radio"/> DR	<input type="radio"/> DAO
<input type="radio"/> Ext 1	<input type="radio"/> I(535nm) + Ext 2	<input type="radio"/> DPM	<input type="radio"/> DNADPH
<input type="radio"/> Ext 2	<input type="radio"/> Ext 1 + Ext 2	<input type="radio"/> DPD	<input type="radio"/> DP700
		<input checked="" type="radio"/> DP515/535	<input type="radio"/> DP515
<b>Analysis Mode</b>		I(535 nm)	P515
<input checked="" type="radio"/> SP-Analysis	<input type="radio"/> Fast Acquisition	Gain	Gain
		<input checked="" type="radio"/> 1 (Low)	<input type="radio"/> 1 (Low)
<input type="checkbox"/> Flux Mode		<input type="radio"/> 5 (High)	<input checked="" type="radio"/> 5 (High)
		Damping	Damping
		<input type="radio"/> 10 $\mu$ s (Low)	<input type="radio"/> 10 $\mu$ s (Low)
		<input checked="" type="radio"/> 1 ms (High)	<input checked="" type="radio"/> 1 ms (High)
		Zero Offset	Zero Offset

## Setup and Connections

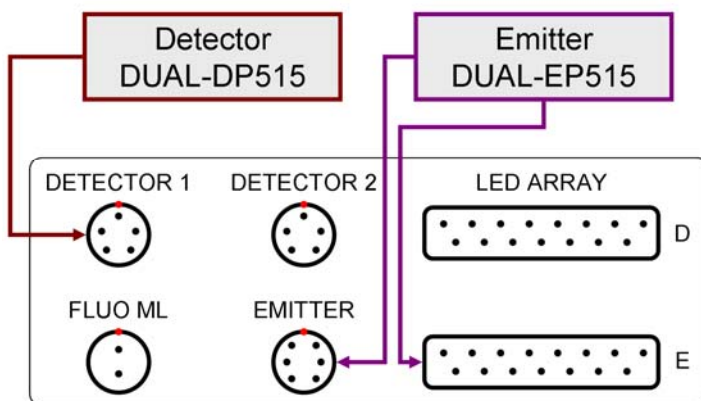
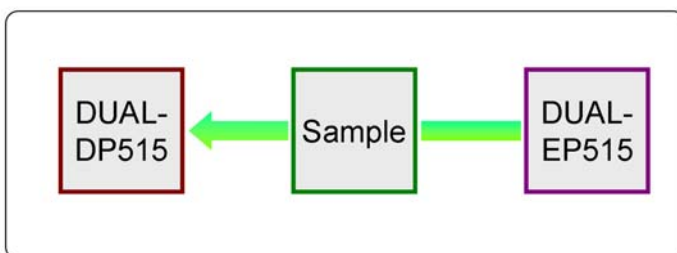


Figure 1. Setup for P515 and I(535) measurements

### 1.3 Setup for parallel assessment of P515 and chlorophyll fluorescence

For concomitant measurements of the P515 change and chlorophyll fluorescence, the DUAL-DP515 detector is connected to the DETECTOR 2 socket of the DUAL-C Power-and-Control-Unit. Then, the fluorescence detector is plugged into socket DETECTOR 1. Chlorophyll fluorescence is detected at an angle of 90 degrees to the optical path of P515 detection (Fig. 2).

When DUAL-DR or DUAL-DB head is used as the fluorescence detector, their cables for modulated fluorescence excitation (fluorescence measuring light) and for actinic light are connected to sockets FLUO ML and LED ARRAY “D”, respectively (Fig. 2). For reasons of electronic configuration, switching on the fluorescence measuring light (F ML) activates not only the fluorescence measuring light in DUAL-DR or DUAL-DB) but also the 535 nm measuring light in the DUAL-EP515.

The MODE-setting window of the DualPAM software must be edited to match the current configuration as shown in Fig. 2. For instance, when a DUAL-DB head is used for fluorescence detection, select DB as detector 1, and DP515 as detector 2 (Fig. 2). Parallel measurements of chlorophyll fluorescence and P515 are enabled by checking “Fluo + P515”, as the DUAL CHANNEL mode.

Make sure that “Fluo Gain” is set to 5 (high). When, after changing of measuring heads, the “Fluo Gain” is set to 1 (low), select 5 (high), point fluorescence emitter away from any objects and carry out a zero offset.

## Mode-setting window of DualPAM software

<b>Measure Mode</b>		<b>Detector Type</b>	
Single Channel	Dual Channel	Detector 1	Detector 2
<input type="radio"/> P515	<input checked="" type="radio"/> Fluo + P515	<input checked="" type="radio"/> DB	<input type="radio"/> N.C.
<input type="radio"/> Fluo	<input type="radio"/> P515 + Ext 1	<input type="radio"/> DR	<input type="radio"/> DAO
<input type="radio"/> Ext 1	<input type="radio"/> Fluo + Ext 2	<input type="radio"/> DPM	<input type="radio"/> DNADPH
<input type="radio"/> Ext 2	<input type="radio"/> Ext 1 + Ext 2	<input type="radio"/> DPD	<input type="radio"/> DP700
		<input type="radio"/> DP515/535	<input checked="" type="radio"/> DP515
<b>Analysis Mode</b>		<b>Fluo Gain</b>	<b>P515 Gain</b>
<input checked="" type="radio"/> SP-Analysis	<input type="radio"/> Fast Acquisition	<input type="radio"/> 1 (Low)	<input type="radio"/> 1 (Low)
<input type="checkbox"/> Flux Mode		<input checked="" type="radio"/> 5 (High)	<input checked="" type="radio"/> 5 (High)
		<b>Damping</b>	<b>Damping</b>
		<input type="radio"/> 10 $\mu$ s (Low)	<input type="radio"/> 10 $\mu$ s (Low)
		<input checked="" type="radio"/> 1 ms (High)	<input checked="" type="radio"/> 1 ms (High)
		<b>Zero Offset</b>	<b>Zero Offset</b>

## Setup and Connections

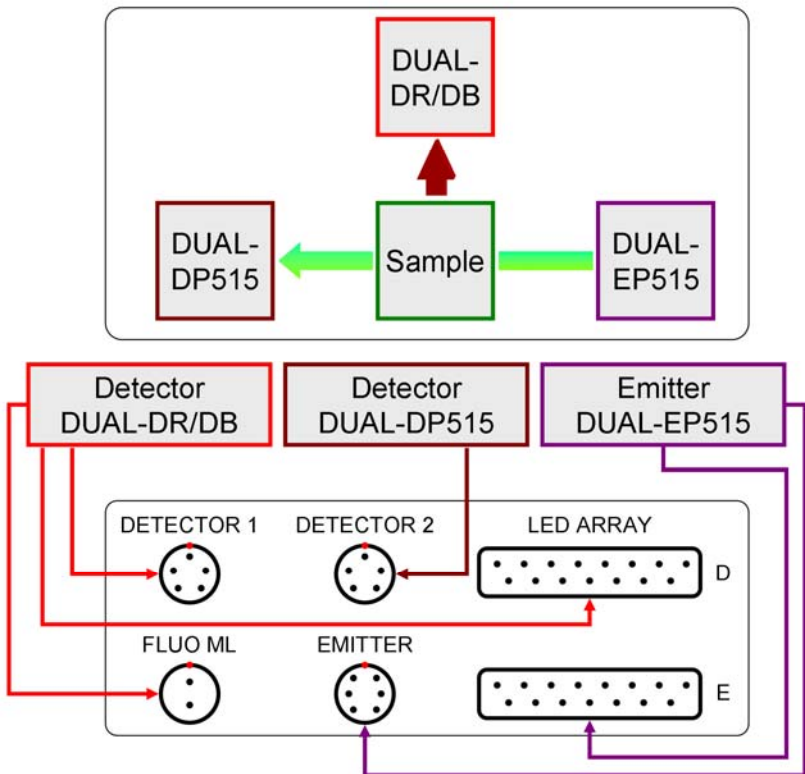


Figure 2. Setup for P515 and chlorophyll fluorescence measurements

## 1.4 Spectra of measuring light, P515 and 535 nm, scattering changes

Figure 3 shows the spectra of light emission by the measuring LEDs of the DUAL-EP515 (upper panel: ML 520, ML 535, and ML 550). In comparison, absorbance difference spectra (treated minus untreated) of the 515 nm electrochromic band shift (also called P515 or P518) and the 535 nm scattering change are shown in the middle and lower panels, respectively. Difference spectra were normalized to the maximum in the green spectral region. Source of data: see overleaf.

Figure 3 represents P515 and 535 nm scattering changes as absorbance data. The P515/535 Emitter-Detector Module, however, does not record absorbance of the sample but the sample's transmission of measuring beams. For example, the P515 signal measured by P515/535 Emitter-Detector Module corresponds to the  $I_{550\text{ nm}}$  minus  $I_{520\text{ nm}}$  signal difference.

For the small changes of P515, however, the variations in  $I_{550\text{ nm}} - I_{520\text{ nm}}$  are positively linearly related to the absorbance changes calculated from the original signals. Therefore, the P515 data measured with the DUAL-PAM-100 can be directly compared with corresponding absorbance data.

Further, the DualPAM software plots the  $I_{535\text{ nm}}$  using an inverted ordinate scale. This representation of the  $I_{535\text{ nm}}$  matches with absorbance data, because the small variations of the 535 nm scattering change *in vivo* are negatively linearly related to the corresponding absorbance changes.

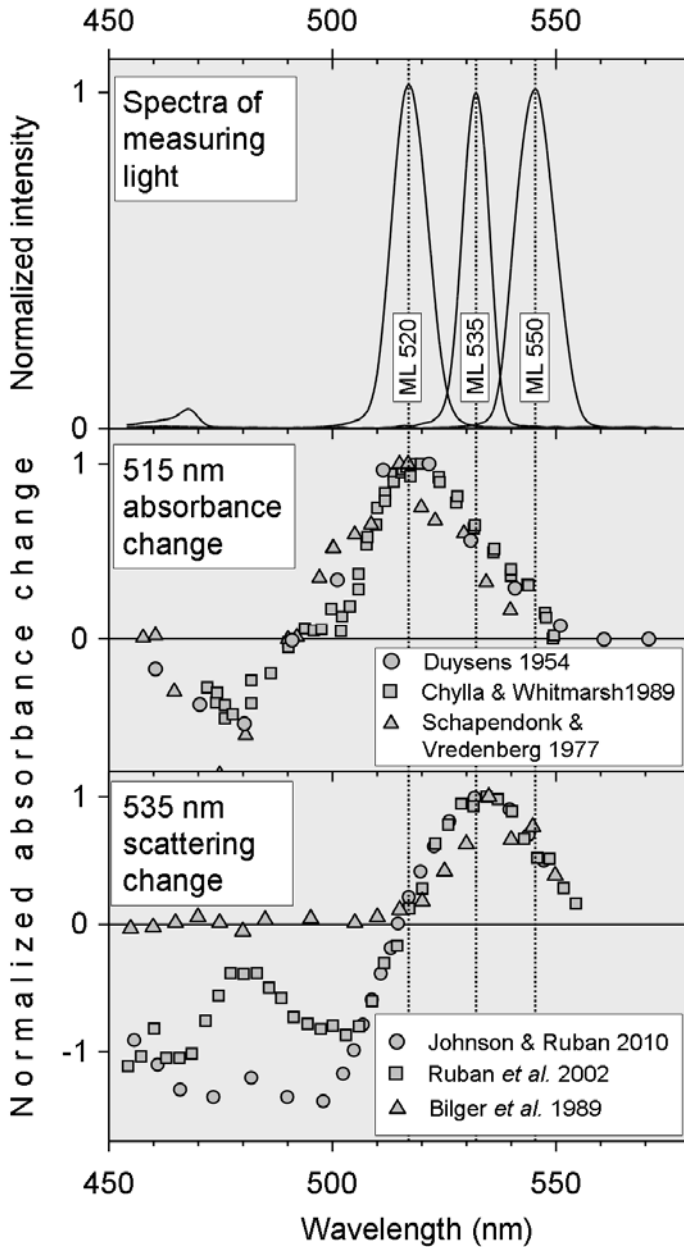


Figure 3. Spectra of measuring light, P515 and 535 nm, scattering changes

### 1.4.1 Data source for Figure 3

Bilger W, Björkman O, Thayer SS (1989) Light-induced spectral absorbance changes in relation to photosynthesis and the epoxidation state of xanthophyll cycle components in cotton leaves. *Plant Physiol* 91: 542-551

Chylla RA, Whitmarsh J (1989) Inactive photosystem II complexes in leaves. Turnover rate and quantitation. *Plant Physiol* 90: 765-772

Duysens LNM (1954) Reversible changes in the absorption spectrum of *Chlorella* upon irradiation. *Science* 120: 353-354

Johnson MP, Ruban AV (2010) Arabidopsis plants lacking PsbS protein possess photoprotective energy dissipation. *Plant J* 61: 283–289

Ruban AV, Pascal AA, Robert B, Horton P (2002) Activation of zeaxanthin is an obligatory event in the regulation of photosynthetic light harvesting. *J Biol Chem* 277: 7785–7789

Schapendonk AHCM, Vredenberg WJ (1977) *Biochim Biophys Acta* 462: 613-621



## 2 Accessory: DUAL-K25 Quartz Glass Cuvette for Suspensions

The DUAL-K25 quartz glass cuvette has been developed to reduce baseline drifts caused by particle settling in suspensions of isolated chloroplasts, unicellular algae and cyanobacteria. Influences of particle settlings are reduced by establishing a vertical optical pathway between emitter and detector units. This way, total light absorption by a suspension stays comparably constant despite particle settling.

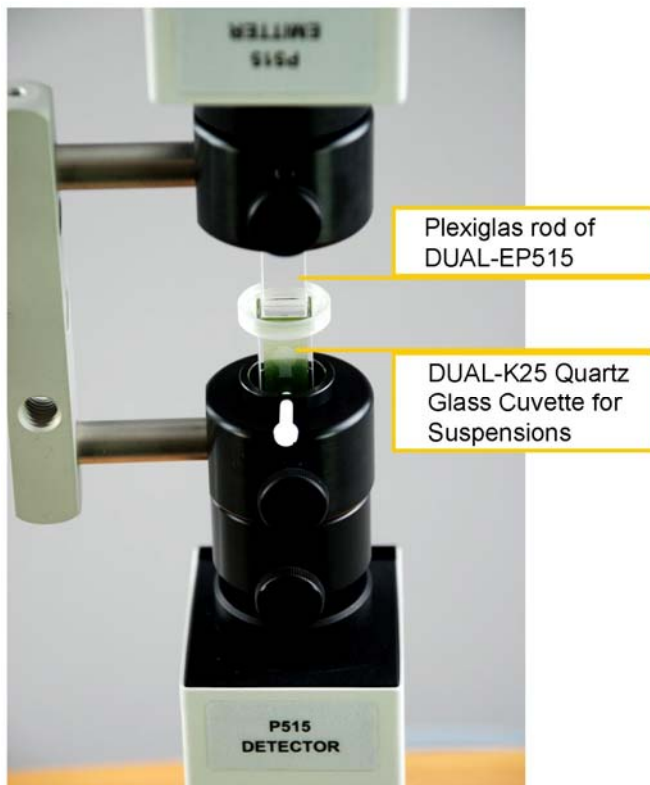


Figure 4. Setup of DUAL-K25 cuvette and P515/535 module



### 3 Index

Artificial Leaf.....	3
Chip-on-board .....	2
DUAL-DP515.....	2
DUAL-K25 cuvette.....	11
Electric field .....	1
Extent of delivery.....	2
Far red light .....	2
Fig. 1. Setup for P515/535 .....	5
Fig. 2. Setup for P515/fluorescence .....	7
Fig. 3. LED emission and in vivo absorbance difference spectra.....	9
Fig. 4. DUAL-K25 cuvette and P515/535 module .....	11
Filter, blue-green .....	2
Multiple turnover flash .....	2
P515/535 Detector Head.....	2
P515/535 Emitter-Detector Module .....	1
P515/P535 Emitter Head.....	2
Setup for P515/535 .....	4
Setup for P515/fluorescence .....	6
Single turnover flash.....	2
Xanthophyll cycle .....	1
Zero offset .....	6



# New accessory for the DUAL-PAM-100: The P515/535 module and examples of its application

Ulrich Schreiber and Christof Klughammer

**Abstract.** The P515/535 module has become available as a new accessory for the DUAL-PAM-100. It serves for simultaneous measurements of the dual beam 550 - 515 nm difference signal, reflecting the electrochromic pigment shift (ECS, P515), and of the single beam 535 nm signal, reflecting "membrane energization" ("scattering" changes). The technical features of the emitter and detector units of the new device are outlined and some typical examples of application are presented. A wide range of frequencies of the pulse modulated measuring light is provided, which can be rapidly altered under software control. In this way, the dark baseline can be assessed (at low frequency) and at the same time rapid induction kinetics can be measured (at high frequency). Due to an exceptionally high signal/noise ratio, the P515 change induced by a single-turnover saturating flash can be reliably measured without averaging. Comparison of the single beam responses at 550, 535 and 515 nm with the dual beam 550 - 515 nm difference signal suggests that the latter is essentially free of "scattering" changes. Besides P515 (ECS), the 550 - 515 nm signal also contains information on the zeaxanthin content, which is reflected in changes of the "dark baseline". Such changes can be reliably evaluated due to the high signal stability, as demonstrated by measurements using an inert sample ("artificial leaf"). Valuable information on the membrane potential ( $\Delta\psi$ ) and proton gradient ( $\Delta\text{pH}$ ) components of the proton motive force (pmf) can be obtained by analyzing light-dark responses of the P515 signal, as previously demonstrated by David Kramer and co-workers (DIRK-analysis). The DUAL-PAM provides a special "P515 Flux" mode of operation, with the actinic light being applied in form of light-dark pulses (1/1 modulation of light/dark). The "P515 Flux" signal provides valuable information on the steady-state rates of coupled electron and proton fluxes *in vivo*.

## Introduction

The DUAL-PAM-100 is the successor of the PAM-100 Chlorophyll Fluorometer introduced in 1985. The "DUAL" stands for simultaneous measurements of two key signals, namely P700 absorbance and chlorophyll (Chl) fluorescence, as indicators of the efficiency of energy conversion in photosystems I and II, respectively. In principle, such measurements have also been possible with the PAM-100, using various accessory components, like special emitter-detector units, single- and multiple-turnover flash lamps, actinic lamps for continuous illumination with blue, red or far-red light and data acquisition system. All these components were miniaturized and optimized for fully computerized operation and integrated in the DUAL-PAM-100 to provide a compact and highly efficient universal photosynthesis measuring system. Besides measuring P700 and Chl fluorescence, this system can be extended by optional emitter-detector modules for assessment of other key photosynthesis parameters, like transthylakoid proton gradient and membrane potential. Generally, the DUAL-PAM-100 is designed for measuring one pulse-modulated **single beam signal** (normally Chl fluorescence excited by red or blue measuring light) and one pulse-modulated **dual beam difference signal** (normally P700 absorbance detected *via* the difference in the transmittance of 870 and 830 nm light).

Here we report on a new P515/535 module, which can serve for *in vivo* assessment of the two components of the proton motive force (pmf), **transmembrane potential**

( $\Delta\psi$ ) and **proton gradient** ( $\Delta\text{pH}$ ), which drive ATP-synthesis *via* the reversible thylakoidal ATP-ase. The membrane potential is measured *via* the "electrochromic pigment absorbance shift" (ECS, P515) characterized by a broad peak at 515 - 525 nm. The  $\Delta\text{pH}$  can be estimated from the light-off kinetics of the P515 signal and independently from the so-called "scattering" signal, characterized by a broad peak at 535 nm.

## Some background on P515.

The P515 signal was discovered by Duysens (1954). Later, Junge and Witt (1968) showed that the P515 signal provides a linear measure of the electric potential across the thylakoid membrane. It is due to an **electrochromic shift** of pigment absorbance, mainly carotenoids and Chl b. For most of the original work on this signal, repetitive flash relaxation spectroscopy has been applied, as high quality signals can be obtained by averaging of the rapid responses, which can be readily distinguished from "non-specific" slow absorbance changes. A saturating single turnover flash induces an extremely rapid signal increase, which reflects the primary **charge separation at PS I and PS II reaction centers**. In intact, dark-adapted samples this rapid rise is followed by a much slower rise phase ( $t_{1/2}$  approx. 20 ms), which is due to the electrogenic step in the so-called **Q-cycle** of the cyt b/f complex (Velthuis 1978), i.e. reflecting secondary inter-system electron transport induced by the primary charge separation. The rate of the ensuing signal decay is highly dependent on the physiological state of the sample, in

particular with respect to membrane integrity and preillumination. The decay is generally accelerated by membrane ion-leakiness and specifically by H<sup>+</sup> conductivity of the reversible thylakoidal ATP-ase. The latter is activated by illumination and linked to transduction of the pmf in ATP-synthesis.

Until recently, relatively little work has been dedicated to the slow changes of P515 absorbance in the minutes time range upon continuous illumination. Klughammer *et al.* (1990) developed a Kinetic LED Array Spectrophotometer that allows to deconvolute a variety of overlapping absorbance changes including P515. Measurements with this system revealed a substantial membrane potential maintained during steady-state illumination, which is highly susceptible to physiological stress. Recently, Sacksteder and Kramer (2000) introduced a special technique for evaluation of the P515 (ECS) signal during steady state illumination: The DIRK (Dark Interval Relaxation Kinetic) technique uses brief (<500 ms) dark intervals to create perturbations in steady-state electron and proton fluxes. The initial rate of the ECS (P515) decay has been suggested to be a linear measure of proton flux rate (efflux *via* ATP-synthase that in the steady-state equals influx *via* protolytic electron transport steps) and consequently of electron flux as well. When longer dark times are given, the ECS (P515) displays complex relaxa-

tion kinetics that have been applied by Cruz *et al.* (2001) for differentiation between  $\Delta\psi$  and  $\Delta pH$  components of the overall pmf (partition analysis). The latter analysis poses considerable demands to the stability and selectivity of the measured signal.

### Some background on the 535 nm signal.

The absorbance change at 535 nm has been attributed to a light induced increase of "light scattering" related to internal acidification of the thylakoids upon  $\Delta pH$  formation (Deamer *et al.* 1966, Heber 1968). It has been used in numerous studies as a convenient semi-quantitative indicator of "membrane energization" and of the  $\Delta pH$  component of the pmf in isolated chloroplasts and intact leaves. It is closely correlated with another semi-quantitative indicator of "membrane energization"/ $\Delta pH$ , namely non-photochemical fluorescence quenching as measured *via* the Chl fluorescence parameters qN (in particular energy-dependent qE) and NPQ (see e.g. Bilger *et al.* 1988). While it has been assumed for a long time that the 535 nm signal is somehow caused by changes in grana stacking, this interpretation recently has been questioned by Ruban *et al.* (2002), who suggest that the increase in 535 nm absorbance is due to a red shift of the zeaxanthin absorbance peak.

## Materials and Methods

### Technical features of the new P515/535 module.

The P515/535 module consists of the emitter unit DUAL-EP515 and the detector unit DUAL-DP515, optical design and outer appearance of which are equivalent to the standard emitter detector units for simultaneous measurements of P700 and Chl fluorescence, DUAL-E and DUAL-DR (or DUAL-DB), respectively. These units can be connected to the same leaf holder or optical unit ED-101US/MD (for measurements with suspensions). Hence, identical optical geometries and light conditions are provided for comparative measurements of the different photosynthesis parameters.

The DUAL-EP515 is connected to the EMITTER socket and the DUAL-DP515 (Fig. 1) to the DETECTOR 1 socket of the DUAL-C Power-and-Control Unit. The P515 signal is measured dual beam as 550 - 515 nm difference signal in analogy to the 870 - 830 nm P700 signal. The 550 nm reference was chosen in order to minimize any contribution of "light scattering" changes to the difference signal. Assuming a Gauss shape absorbance peak at 535 nm, the absorbance changes at 550 and 515 nm should be about equal (see related data below under Results). On the other hand, in contrast to 550 nm at 515 nm there is substantial absorbance by zeaxanthin (broad absorbance band), so that the 550-515 nm signal not only reflects P515, but zeaxanthin changes as well. These two signal components, however, can be readily distinguished by the fact that zeaxanthin changes are extremely slow. The dark-light and light-dark induced



Fig. 1. Emitter unit DUAL-EP515 and detector unit DUAL-DP515

550 – 515 nm changes in the sub-s to min time range are almost exclusively due to P515, whereas the slower changes ranging from minutes to hours reflect zeaxanthin formation (*via* violaxanthin and antheraxanthin deoxidation) and zeaxanthin epoxidation, respectively.

**Specifications of the DUAL-EP515.** The DUAL-EP515 features an array of 8 white LEDs equipped with interference filters. While the 550 nm measuring light is derived from 3 LEDs with 3 individual 550 nm interference filters (5.5 nm HBW), 4 LEDs equipped with 4 individual 515 nm interference filters (8.5 nm HBW) provides 520 nm measuring light which was actually used to measure the “515 nm” absorption changes. A single LED with a 535 nm interference filter (5.5 nm HBW) gives 535 nm measuring light. These LEDs, which are arranged in a ring, are focused *via* a central 6.5 mm hole in a chip-on-board (COB) LED array (featuring 635 nm Power-LEDs) on a 10 x 10 mm Perspex rod, which serves for mixing the various light qualities and guiding the randomized light to the sample. In addition, a single 740 nm LED equipped with a 1 mm RG9 filter in the center of the LED array serves for far-red illumination with preferential excitation of PSI. The COB array consists of 24 Power-LEDs, which for short time intervals can be driven with high currents (up to 1.5 A). It provides continuous actinic illumination as well as Multiple Turnover and Single Turnover flashes.

**Specifications of the DUAL-DP515.** The DUAL-DP515 features a 10x10 mm PIN-photodiode equipped with a 2 mm blue-green filter (BG39) transmitting the 520, 535 and 550 nm measuring light but blocking the 635 nm actinic and 730 nm far-red light. A 10 x 10 mm Perspex rod guides the transmitted measuring light from the sample to the BG39 filter in front of the PIN-photodiode, which is mounted on a printed circuit board, where the pulse-modulated signals are preamplified. It should be noted that for the sake of maximal signal/noise the DUAL-DP515, in contrast to the DUAL-DR (or DUAL-DB), is not equipped with a COB LED array for actinic illumination. Light-induced responses in P700/Chl fluorescence can be compared at identical intensity settings with the corresponding responses in P515/“Scattering” when the LED Array cable of the DUAL-DR (or DUAL-DB) is disconnected.

**Calibration of the P515 signal.** The DUAL-PAM-100 applies a special modulation technique for high-sensitivity measurements of dual beam difference signals. The DualPAM software supports automated balancing of the two transmitted beams *via* adjusting the LED currents (“Balance”, Fig. 2), so that the difference signal is close to zero. The resulting LED current ratio, which depends on the spectral properties of the sample, is displayed. After “Balance” the changes in the difference signal can be detected with high sensitivity without any risk of amplifier saturation. The single signals, which may be more than 1000x larger than the recorded signal changes, can be determined with the help of a special calibration rou-

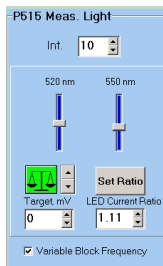


Fig. 2. Balancing the 520 and 550 nm beams - screenshot.

time, which involves a defined transient decrease of the 515 nm signal with respect to the 550 nm signal (*via* corresponding decrease in LED current). The original difference signals are measured in Volt units, which are transformed into  $\Delta I / I$  units with the help of the calibration routine. The “I” corresponds to the amplitude of the balanced overall signals. In Fig. 3, an example of calibration is shown for a tobacco leaf at Measuring Light intensity setting 10. In this example the  $\Delta I / I$  corresponding to the depicted bar amounts to  $3.57 \times 10^{-3}$  units and an overall signal  $I = 26.8 \text{ V}$  is calculated.

In principle, it is also possible to measure the single beam 515 nm and 550 nm signals (see example below under Results), when the LED current slide of the other (normally compensating) LED is moved to zero. In this case, Gain 1 (Low Gain) has to be applied to avoid amplifier saturation.

**Recording of the single beam 535 nm signal.** The single beam 535 nm signal is recorded in form of a non-compensated signal in analogy to Chl fluorescence. Just like the “Fluo” signal, it saturates at about 5.5 V. Therefore, in contrast to the dual beam 550 - 515 nm signal, it is recorded at Low Gain (Fig. 4A), so that a high setting of Measuring Light (ML) intensity can be applied, for the

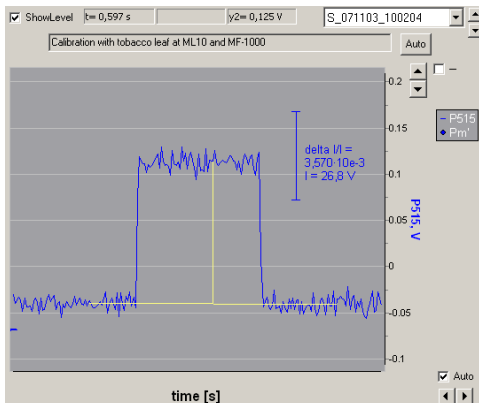


Fig. 3. Definition of the signal scale in  $\Delta I / I$  units for a tobacco leaf using the automated calibration routine.

sake of a high signal/noise ratio. As "scattering" changes are relatively slow, High Damping can be applied, thus further improving the signal/noise. For display of kinetic recordings the 535 nm signal is inverted, so that a "scattering" increase, which corresponds to a signal decrease, is displayed as an upwards transient, in analogy to the P515 signal.

### Means to prevent preillumination effects of the Measuring Light.

In contrast to measurements of P700 *via* the 870-830 nm signal, where near-infrared ML can be applied at very high intensity without any actinic effects, measurements of P515 and "scattering" are complicated by a substantial fraction of the green measuring light being absorbed by photosynthetic pigments. The DUAL-PAM software provides a number of means for preventing significant preillumination effects of the ML without sacrificing signal quality. This particularly applies to the rapid P515 change induced by a single turnover saturating flash (ST flash), the kinetics of which are extremely sensitive to preillumination (due to light activation of the reversible ATPase). In this type of measurements advantage is taken of the fact that the LED driven ML can be rapidly switched on. When "Auto ML on" is active (Fig. 4B), upon Start of a Fast Kinetics recording the ML is automatically switched on, so that before the measurement the sample can be kept in absolute darkness. Furthermore, use of the so-called "Variable Block Frequency" can be made (Fig. 4C). ML pulses are applied in the form of 30  $\mu$ s "pulse

blocks" (containing 12 pulses) separated by more or less long "dark blocks". The signal level may be monitored with negligibly small actinic effect at low "Block Frequency" (MF-low ranging from 1-1000 Hz), e.g. for signal calibration or for assessment of the absolute level. The actual measurement of light-induced changes can be carried out at high pulse frequency (MF-high ranging from 100-10000 Hz), depending on the required time resolution. For maximal time resolution the ML pulse frequency can be switched to MF-max, corresponding to about 400 kHz, with the help of pre-programmed Fast Trigger files, which define triggering of the various light sources and of MF-max with 2.5  $\mu$ s resolution.

### Signal/noise ratio and signal stability.

For reliable *in vivo* measurements with intact leaves particularly high demands on signal/noise ratio and signal stability are made. The signal amplitude depends on the intensity of the measuring light (settings 1 - 20), the amplifier gain (settings 1 and 5) and the transmittance of the sample. While with leaf samples and isolated chloroplasts the noise is mainly determined by the electronics, in the case of algae suspensions the "sample noise" (movement of cells or groups of cells) contributes substantially to overall noise. Signal stability, which is particularly important for long term recordings, is affected by instrumental and sample parameters as well. In particular, signal drift can be caused by temperature related intensity changes of the pulse-modulated LED measuring light and by optical changes of the sample (e.g. caused by fluctuations of the leaf water status or settling of algae). As long as these changes apply equally to the 515 nm and 550 nm signals, they should be eliminated in the 550 - 515 nm difference signal.

The noise amplitude of a given recording is a function of the number of "Point Averages", the Measuring Light Frequency (MF-low, MF-high and MF-max) and of the Damping (time constants 10  $\mu$ s for Fast Kinetics or 1 ms for Slow Kinetics). The DUAL-PAM-100 software offers an exceptionally wide range of data acquisition rates. Fast Kinetics data can be acquired at rates ranging from 2.5  $\mu$ s to 1000  $\mu$ s/point, with up to 128000 points being saved. Slow Kinetics can be recorded at 1 - 100 ms/point with up to 1024000 points (Fig. 4D); maximal recording time more than 28 hours). In practice, kinetics normally are measured with acquisition rates much higher than required for time resolution of the recorded kinetics and at a later stage, upon viewing and analyzing the data, suitable "Point Averaging" is applied for noise reduction (Fig. 4E). For example, when the flash-induced P515 change (consisting of an instantaneous rise followed by slower changes in the ms time range) is measured at 2.5  $\mu$ s/point and MF-10000, no kinetic information is lost, if 100 "Point Averages" are applied (see example under Results in Fig. 6). In this particular example, the noise amplitude amounted to  $10^{-4}$   $\Delta I/I$  units with 100 "Point Averages", as compared to  $3 \times 10^{-3}$   $\Delta I/I$  units without "Point Averaging". In view of a flash-induced signal change in

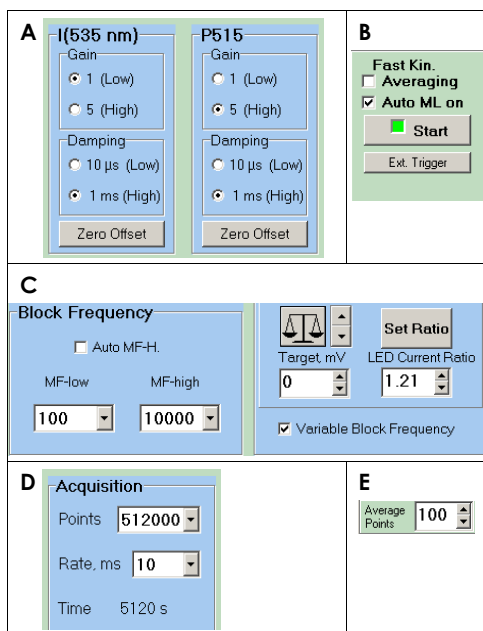


Fig. 4. Parameter adjustment: screenshots.



the order of  $15 \times 10^{-3} \Delta I/I$  units, this noise may be considered rather low.

A non-ambiguous assessment of the instrument related noise and signal "drift" is possible using an "artificial leaf". For this purpose, a plastic filter sheet with a transmittance spectrum in the green region similar to a green leaf (Roscolux #01: Light Bastard Amber) was sandwiched between two layers of highly scattering transparent paper. With this "artificial leaf #01" the signal amplitudes at 515, 535 and 550 nm were very similar to those of a normal tobacco leaf.

Figure 5 shows the 550 - 515 nm signal measured with the "artificial leaf #01" at increasing Measuring Light Frequencies (MF). At higher MF settings on one hand the noise is reduced and on the other hand a negative signal drift is induced. While the noise at MF20 amounts to about  $3 \times 10^{-3} \Delta I/I$  units, it decreases to about  $10^{-4} \Delta I/I$  units at MF1000 and  $3 \times 10^{-5} \Delta I/I$  units at MF10000. At MF10000 and ML Intensity 10 the integrated photosynthetically active radiation (PAR) of the ML amounts to  $4 \mu\text{mol} \mu\text{mol quanta} \cdot \text{m}^{-2} \cdot \text{s}^{-1}$ . In most applications, the actinic effect of the ML may be considered negligibly small, when applied below MF5000. Up to MF1000 the signal drift over 5 min time period does not exceed  $10^{-4}$

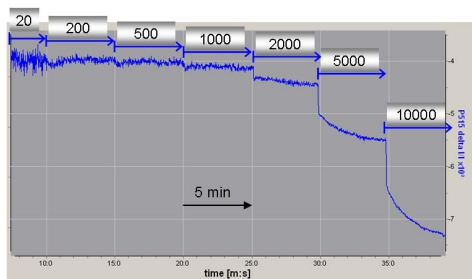


Fig. 5. Effect of increasing Measuring Light Frequency (MF) on the noise and stability of the 550-515 nm signal using the "artificial leaf #01". ML intensity 10. Points averaged: 30.

$\Delta I/I$  units. Upon switching to higher MF settings there is a biphasic decline of the signal, which at MF2000 amounts to about  $3 \times 10^{-4}$  and at MF10000 to about  $2 \times 10^{-3} \Delta I/I$  units over 5 min time periods.

High MF settings normally are applied for the recording of sub-s Fast Kinetics only, where slow signal drifts are of no concern and where the actinic effect of the ML can be ignored.

## Results and Discussion

### Rapid P515 relaxation kinetics induced by single-turnover saturating flash

In Fig. 6 typical recordings of P515 absorbance changes induced by a saturating single turnover flash in an intact tobacco leaf are presented. Fig. 6A shows the response after 1 h dark adaptation. A rapid rise of  $14 \times 10^{-3} \Delta I / I$  units is followed by a distinct slow rise phase, with the overall response peaking at about 20 ms, before a slow decline sets in. In Fig. 4B the response of the same sam-

ple is shown after 10 min illumination at  $100 \mu\text{mol quanta} \cdot \text{m}^{-2} \cdot \text{s}^{-1}$  followed by 4 min dark adaptation. The preillumination has dramatically speeded up the decay of the flash induced P515 change and the slow rise phase has disappeared. The rapid rise phase is diminished by about 13%.

It may be emphasized that the P515/535 module allows to assess the essential parameters of the flash induced P515 change with single recordings, i.e. there is no need for


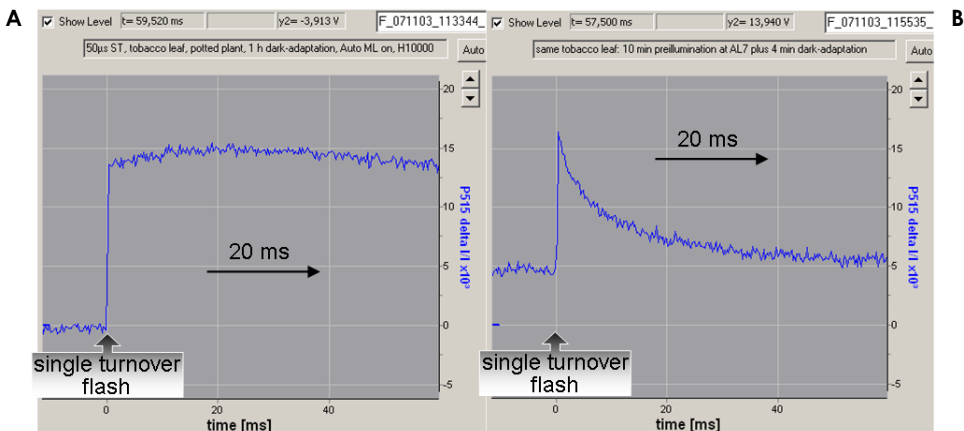


Fig. 6. Typical recordings of P515 changes induced by saturating single turnover flashes in tobacco leaf. Single recordings of Fast Kinetics (no curve averaging) using Auto ML on and MF-10000. Display with 100 points averaging. A) Dark-adapted for 1 h. B) Preilluminated for 10 min at  $100 \mu\text{mol quanta} \cdot \text{m}^{-2} \cdot \text{s}^{-1}$  followed by 4 min dark.

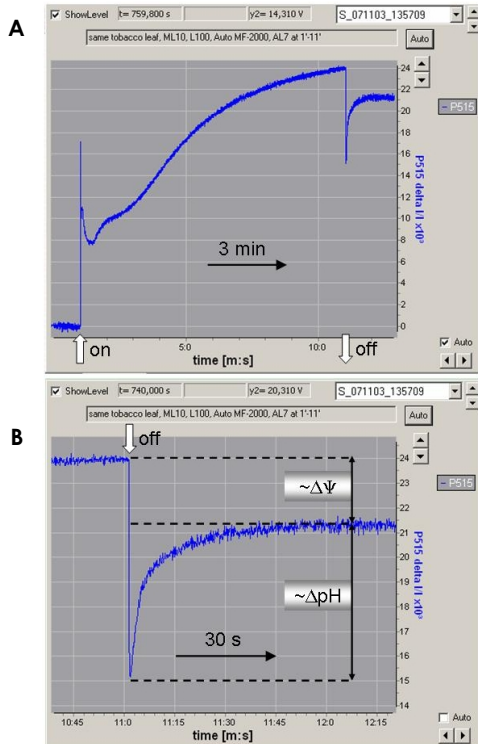

**Fig. 7.** Slow P515 changes induced by actinic illumination (10 min at  $100 \mu\text{mol quanta/m}^2\text{s}$ ) and return to darkness. Same tobacco leaf as in experiment of Fig. 4 after dark-adaptation. Use of Auto MF-High, switching from MF 200 to 2000. A) Complete recording of light-on and light-off responses. B) Enlarged display of light-off response with indication of estimated  $\Delta\psi$  and  $\Delta\text{pH}$  components of pmf.

curve averaging. This is particularly important for measurements of the slow rise phase, which is affected by the preillumination caused by repetitive measurements. Preillumination activates the reversible ATP-ase in the thylakoid membrane, thus increasing the  $\text{H}^+$  conductivity of the membrane. The decay of the P515 signal reflects the relaxation of the flash induced electric field (created by charge separation in the two photosystems and electrogenic electron transport the Q-cycle at the cyt b/f complex) by  $\text{H}^+$  efflux *via* the  $\text{H}^+$  channel of the ATP-ase. A functionally intact photosynthetic apparatus is characterized by a slow decay after dark-adaptation (high membrane integrity) and a fast decay after illumination (high ATP-ase activity).


While in principle besides the electrochromic shift (P515, ECS) other absorbance changes (like zeaxanthin formation or "scattering") can contribute to the 550 - 515 nm difference signal, these alternative changes are much slower and, hence, the rapid flash-induced changes (in

the sub-s time range) may be considered to reflect P515 only. This is not true for slow changes in the time range of min (see below).

### Slow dark-light-dark induction transients of the 550-515 nm signal

Fig. 7A shows the recording of dark-light and light-dark induced slow P515 changes in an intact tobacco leaf. Before this measurement, the leaf was kept for several hours in darkness, resulting in low zeaxanthin content. Under these conditions, the light-induced signal increase not only reflects an increase of the membrane potential (ECS), but also formation of zeaxanthin. The latter is a slow process, which sets in upon accumulation of protons within the thylakoid lumen. The relative extent of zeaxanthin formation can be judged from the increase of the "dark baseline" apparent after light-off. Due to the high stability of the 550 - 515 nm difference signal, the observed slow changes in the "dark baseline" can be reliably evaluated in terms of **reversible changes of zeaxanthin content**. After the light induced increase of the "dark baseline" shown in Fig. 7A, there was a very slow signal decline in the dark ( $t_{1/2}$  approximately 20 min) reflecting zeaxanthin epoxidation (not shown in the figures).

The **rapid light-off response**, which is depicted in more detail in Fig. 7B, reflects  $\text{H}^+$  efflux from the lumen to the stroma *via* the thylakoid ATP-ase. The rapid signal decline is followed by a biphasic signal increase to an apparent "dark baseline". According to David Kramer and co-workers (Kramer and Sacksteder 1998, Cruz *et al.* 2001) the **relative amplitudes of  $\Delta\psi$  and  $\Delta\text{pH}$**  can be estimated from the characteristic levels observed during the light-off response, as indicated in Fig. 7B. The difference between the steady state signal and the "dark baseline" reflects a substantial  $\Delta\psi$  during steady state illumination. The "undershoot" below the "dark baseline" is considered a measure for the steady state  $\Delta\text{pH}$ . During prolonged illumination, proton accumulation within the lumen is

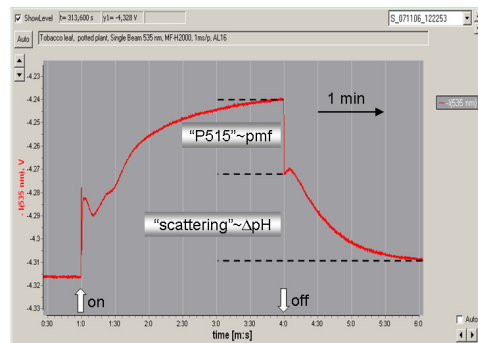

**Fig. 8.** Single beam 535 nm recording of dark-light-dark induction transients. Tobacco leaf of potted plant illuminated for 3 min at  $800 \mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$ . The relative extents of the total proton motive force (pmf) and the  $\Delta\text{pH}$  component are estimated from the rapid and slow light-off responses, respectively (see text).

slowly compensated by anion influx. When upon light-off the accumulated protons are rapidly released from the lumen to the stroma *via* the ATP-ase, there is a sudden excess of negative charges at the internal side of the membrane, resulting in an inversed ECS. The inversed ECS slowly relaxes again during the following dark period, as anion concentrations of lumen and stroma equilibrate.

### Dark-light-dark induction transients of the 535 nm signal

The single beam 535 nm signal shows complex dark-light-dark induction transients, which at low actinic intensities are dominated by the electrochromic shift and at high light intensities by the "scattering" change. Fig. 8 shows a typical recording with a tobacco leaf illuminated for 3 min at  $800 \mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$ . While the rapid on/off transients are due to the ECS (P515), the slow rise and decay phases reflect "scattering" changes, paralleling formation and decline of the  $\Delta\text{pH}$  (Bilger *et al.* 1989).

It should be noted that the amplitude of the rapid light-off

response of the 535 nm signal cannot serve as a measure of the transmembrane  $\Delta\psi$ , as has been assumed in numerous previous studies. In contrast to the rapid light-off response of the 550-515 nm signal described above (see Fig. 7b), it is not possible to measure the "dark baseline", as the slow ECS signal increase due to anion efflux is overlapped by the slow signal upon relaxation of the "scattering" signal.

### Dual beam compared to single beam responses

The main purpose of measuring P515 dual beam as 550-515 nm difference signal is the elimination of the overlapping "scattering" change. The following data suggest that "scattering" changes indeed are successfully suppressed. In Fig. 9 the single beam responses at 515, 550 and 535 nm are compared with the dual beam response (550-515 nm). Single beam 515 and 550 nm signals can be recorded, when after balance of the difference signal one of the two beams is turned off (by setting the LED current slide to zero).



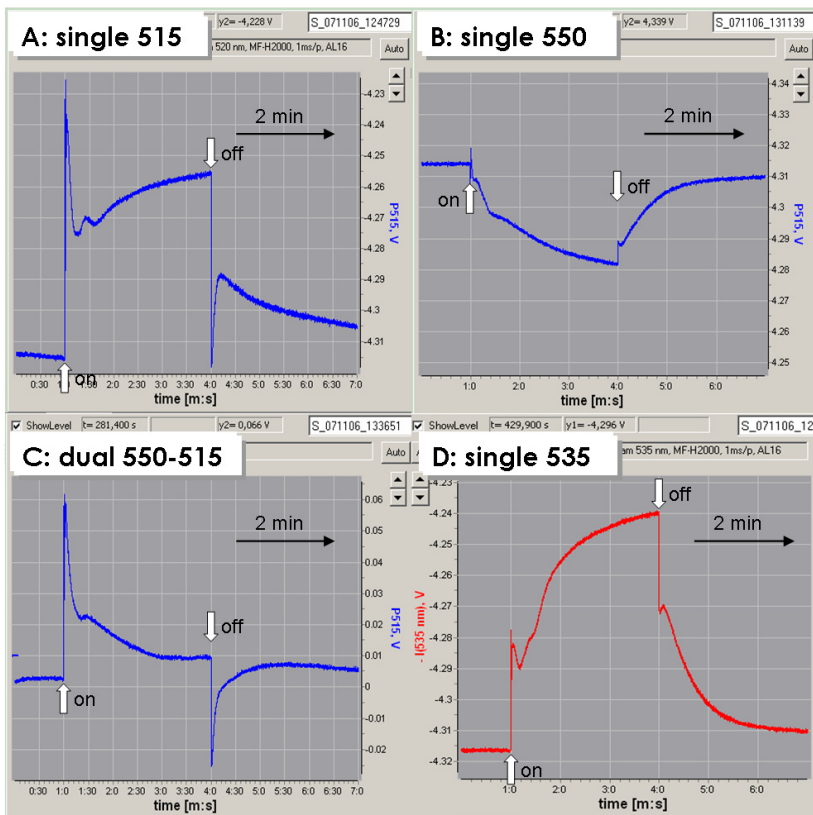


Fig. 9. Comparison of single beam and dual beam responses. Tobacco leaf of potted plant illuminated at  $800 \mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$ . 8 min dark adaptation between consecutive measurements.

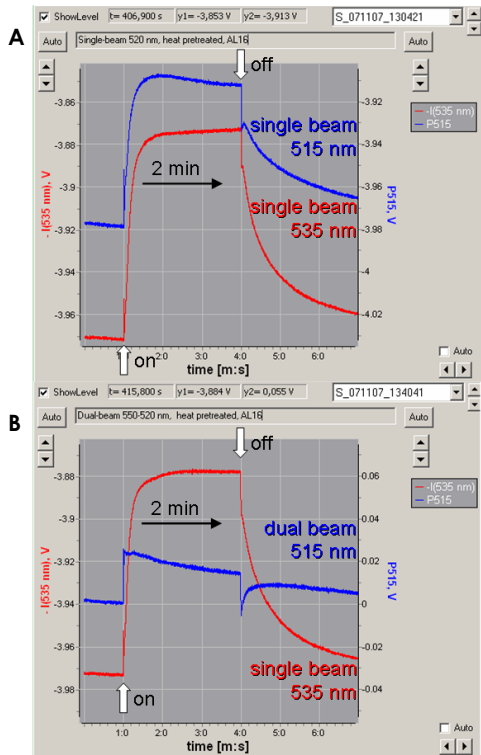

**Fig. 10.** Simultaneous recording of dark-light-dark transients of single beam and dual beam signals of tobacco leaf pretreated for 5 min at 41°C. 30 min recovery at room temperature (20°C). A) Comparison of single beam 515 and 535 nm signals. B) Comparison of dual-beam 550-515 nm with single beam 535 nm signal.

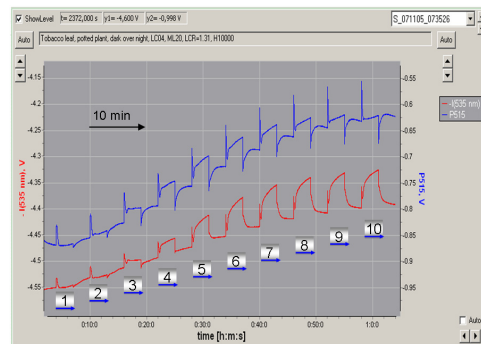

**Fig. 11.** Simultaneous recording of P515 and 535 nm signals of a dark-adapted tobacco leaf during the course of 3 min illumination periods at increasing actinic intensities separated by 3 min dark periods. Actinic intensity sequentially was increased in ten steps from 10 to 1600  $\mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$ . Plant dark-adapted over night before start of Script-file recording.

While the single beam 515 and 550 nm responses display essentially the same slow "scattering" changes that dominate the 535 nm response, these are not apparent in the 550-515 nm difference signal.

Relatively mild heat pretreatment is known to substantially stimulate light induced "scattering" changes and to suppress the light induced P515 change (Bilger and Schreiber 1990). Fig. 10A shows the simultaneously measured single beam responses at 515 and 535 nm of a tobacco leaf pretreated for 5 min at 41°C. Compared with the corresponding traces in Fig. 9, the scattering response is considerably speeded up and its amplitude doubled. At the same time the rapid on/off P515 responses are lowered by about a factor of 3. As shown in Fig. 10B, even under these highly unfavorable conditions, no scattering changes are apparent in the 550 - 515 nm difference signal.


### Long term P515 and 535 nm recording for assessment of light response characteristics

After stabilization of the Measuring Light LEDs (see section on signal stability under Materials and Methods), the P515 and 535 nm signals are sufficiently stable for assessment of long term signal changes in the hour time range. Prerequisite for such measurements is a stable water status of leaf samples. In Fig. 11 an example of a lengthy recording of the light-intensity dependence of the dark-light-dark responses is presented. The DualPAM software allows to program so-called Script-files, consisting of a list of commands, which are sequentially carried out at defined times by the instrument. In principle, all operations that can be carried out manually can also be programmed in a Script-file. In this way, even quite extended and complex experiments can be reliably reproduced.

While a detailed discussion of the complex information contained in the data of Fig. 11 would be out of scope of the present report, the following aspects shall be pointed out:

- Both signals show a substantial rise of the "dark baseline", which may be assumed to reflect formation of zeaxanthin. After light-off, the "dark baseline dark" slowly returns to its original level with  $t_{1/2} = 20$  min (not in the figures).
- Surprisingly zeaxanthin formation already sets in at the lowest actinic intensity (step 1,  $10 \mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$ ), i.e. under clearly light-limited conditions.
- At low actinic intensity P515 and 535 nm signals display very similar changes; both signals are dominated by the ECS, while no scattering change is apparent.
- The steady state  $\Delta\psi$ , as determined by the approach depicted in Fig. 7B, is negligibly small up to about  $60 \mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$  (illumination step 3), then increases up to  $300\text{-}400 \mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$  (illumination steps 6-7), from whereon it declines again to zero.

e) In contrast, the  $\Delta pH$  component of the pmf continuously increases up to the maximal applied intensity. This finding is in agreement with the notion that the  $\Delta pH$  (in contrast to the  $\Delta \psi$ ) plays a key role in the protection against damage by excess light.

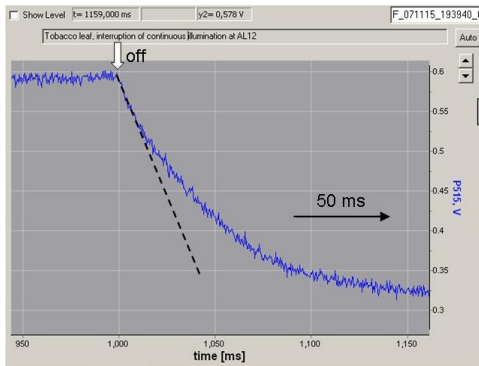

**Fig. 12.** Light-off response of the P515 signal in the steady-state with continuous illumination at  $350 \mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$ . Average of 4 recordings involving 2 s dark-periods separated by 30 s illumination. Tobacco leaf on potted plant.


**Estimation of coupled *in vivo*  $e^-$  and  $H^+$  transport rates via "P515 flux" measurements**

David Kramer and co-workers have presented extended work on the application of "Dark Induced Relaxation Kinetics" (DIRK analysis) for assessment of the pmf and its relative partitioning into  $\Delta \psi$  and  $\Delta pH$  (Kramer and Sacksteder 1998, Kanazawa and Kramer 2002, Cruz *et al.* 2004). This work has confirmed and extended the results of basic research by numerous other researchers (e.g.

Witt, Junge, Rumberg, Joliot, Delosme, Bouges-Bocquet, Velthuys, Vredenberg, van Kooten, Hope, Cramer, Crofts, Whitmarsh, Rich, Garab, Hind, Crowther, Shahak, Morita) showing that the P515 signal carries fundamental information on the rates of coupled electron and proton transport. During continuous illumination a flux equilibrium is established between formation of  $\Delta \psi$  *via* charge separation at the reaction centers and relaxation of  $\Delta \psi$  caused by  $H^+$  efflux *via* the ATP-ase. When illumination is interrupted, formation of  $\Delta \psi$  instantaneously is stopped (except for the dark electrogenic Q-cycle), whereas  $H^+$  efflux first continues at unchanged rate. The initial rate of the P515 light-off response may be taken as a measure of the  $H^+$  efflux rate during illumination, which is a function of the pmf (driving force) and the proton conductivity of the ATP-ase.

Fig. 12 shows an example of measuring the light-off response of the P515 signal in a tobacco leaf with the DUAL-PAM-100. The system was programmed to interrupt continuous illumination at  $350 \mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$  for 2 s and to measure the Fast Kinetics with 4 averages. The initial maximal rate is observed up to about 20 ms after light-off.

The DUAL-PAM-100 software offers the so-called "Flux Mode" for P515 and P700 measurements, which allows continuous assessment of the initial rates of P515 relaxation and P700 re-reduction, respectively. When the "Flux Mode" is active, the actinic light is applied in form of light-dark pulses (1/1 modulation of light/dark). The modulation rate is determined by the Slow Kinetics Acquisition Rate. If e.g. an Acquisition Rate of 20 ms/p is selected, 20 ms AL-pulses are alternating with 20 ms dark-pulses. During steady-state illumination, the increase of the P515 signal (or  $\Delta \psi$ ) during an AL-pulse

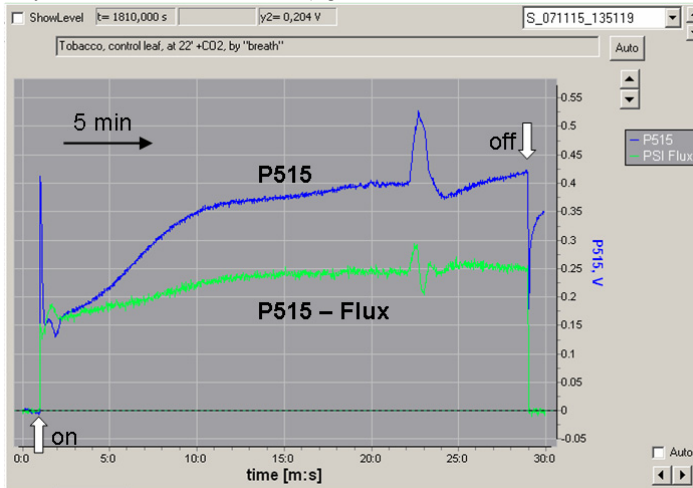

**Fig. 13.** Dark-light induction curves of P515 (blue) and P515 Flux (green) signals. Slow Kinetics Acquisition Time 20 ms/p. Dark-adapted tobacco leaf (attached, potted plant) in standard leaf-holder (without control of gas composition) illuminated at a mean actinic intensity of  $300 \mu\text{mol quanta}\cdot\text{m}^{-2}\cdot\text{s}^{-1}$ . At time 22 min the  $CO_2$  supply was temporarily increased by blowing breath in the direction of the leaf from 1 m distance.

equals the decrease of the P515 signal (or  $\Delta\psi$ ) during a dark-pulse. The repetitive signal changes are processed by the DUAL-PAM-100 to give a continuous "P515 Flux" signal, which can be displayed in parallel with the original P515 signal.

In Fig. 13 a "P515 Flux" measurement of a dark-light induction curve with 20 ms/p Acquisition Time is presented. While phenomenologically the flux signal displays some similarity to the original P515 signal, distinguishing features are:

- a) In contrast to the P515 signal, the flux signal is zero in the absence of actinic illumination. According to its definition, the measured flux corresponds to the charge separation rate in PS I and PS II, which is zero in the dark.
- b) As the flux signal is derived from the rapid AL-off/on changes, it does not contain any contribution of the slow zeaxanthin (and scattering) changes.

When the CO<sub>2</sub> concentration is transiently increased, this leads to a pronounced transient stimulation of the P515 signal. Simultaneously, a damped oscillation of the flux

signal is induced, with initial stimulation being followed by suppression etc. It appears likely that the suppression is due to a transient depletion of ADP, which temporarily slows down the H<sup>+</sup> efflux rate *via* the ATP-ase. The observed oscillation is expression of a disturbance of the flux equilibrium between H<sup>+</sup> influx (from stroma to lumen side of the thylakoid membrane) coupled to photosynthetic electron flow and H<sup>+</sup> efflux *via* the ATP-ase. The oscillation and the eventually reestablished stable flux rate are the result of complex regulatory processes, leading to equilibration of the rates of primary charge separation (with consequent e<sup>-</sup> transport and H<sup>+</sup> influx) and secondary H<sup>+</sup> efflux.

In principle, it is possible to calibrate the "P515 Flux" signal in absolute units of H<sup>+</sup>/(mg Chl-time) or e<sup>-</sup>/(mg Chl-time), based on the amplitude of the rapid phase of the P515 change induced by a single turnover saturating flash. The latter corresponds to the transmembrane field created by the movement of two electrons from inside to outside (one electron each in PS I and PS II). Details on the calibration procedure as well as comparative measurements of the "P700 Flux" signal will be presented in a forthcoming communication.

## References

- Bilger W and Schreiber U (1990) Chlorophyll luminescence as indicator of stress-induced damage to the photosynthetic apparatus. Effects of heat-stress in isolated chloroplasts. *Photosynth Res* 25: 161-171
- Bilger W, Björkman O, Thayer SS (1989) Light-induced spectral absorbance changes in relation to photosynthesis and the epoxidation state of xanthophyll cycle components in cotton leaves. *Plant Physiol* 91: 542-551
- Bilger W, Heber U and Schreiber U (1988) Kinetic relationship between energy-dependent fluorescence quenching, light scattering, chlorophyll luminescence and proton pumping in intact leaves. *Z Naturforsch* 43c: 877-887
- Cruz JA, Avenson TJ, Kanazawa A, Takizawa K, Edwards GE and Kramer DA (2004) Plasticity in light reactions of photosynthesis for energy production and photoprotection. *J Exp Botany* 56: 395-406
- Deamer DW, Crofts AR and Packer L (1966) Mechanisms of light-induced structural changes in chloroplasts. I. Light-scattering increments and ultrastructural changes mediated by proton transport. *Biochim Biophys Acta* 131: 81-96
- Duysens LNM (1954) Reversible changes in the absorption spectrum of *Chlorella* upon irradiation. *Science* 120: 353-354
- Heber U (1969) Conformational changes of chloroplasts induced by illumination of leaves in vivo. *Biochim Biophys Acta* 180: 302-319
- Junge W and Witt HT (1968) On the ion transport system in photosynthesis: investigations on a molecular level. *Z Naturforsch* 23b: 244-254
- Kanazawa A and Kramer DA (2002) In vivo modulation of non-photochemical quenching (NPQ) by regulation of the chloroplast ATP synthase. *Proc Natl Acad Sci USA* 99: 12798-12794
- Klughammer C, Kolbowski J and Schreiber U (1990) LED array spectrophotometer for measurement of time resolved difference spectra in the 530–600 nm wavelength region. *Photosynth Res* 25: 317-327
- Kramer DA and Sacksteder CA (1998) a diffused-optics flash kinetic spectrophotometer (DOFS) for measurements of absorbance changes in intact plants in the steady-state. *Photosynth Res* 56: 103-112
- Ruban AV, Pascal AA, Robert B, Horton P (2002) Activation of zeaxanthin is an obligatory event in the regulation of photosynthetic light harvesting. *J Biol Chem* 277: 7785–7789
- Sacksteder CA and Kramer DA (2000) Dark interval relaxation kinetics of absorbance changes as a quantitative probe of steady-state electron transfer. *Photosynth Res* 66: 145-158
- Velthuys BR (1978) A third site of proton translocation in green plant photosynthetic electron transport. *Proc Natl Acad Sci USA* 76: 2765-2769