

Mounting of Battery Holder #2

DIVING-PAM-II Accessory.

Walz Order Code 000246004613, for 134 x 69 x 36.5 mm,
8.0 V 3.0 Ah, Multipower MP 3-8 battery

Contents

1. Safety Instructions	1
2. Extent of Delivery	2
3. Removing Battery Holder	2
4. Mounting Battery Holder #2	3
5. Connecting the Battery	5

1. Safety Instructions

Do not deviate from the instructions given below and obey the sequence of manipulations.

Avoid short-circuiting battery poles or short-circuiting battery and housing. Watch open ends of adapter cables.

Never connect red and black electric components (wires or battery poles).

When disconnecting the battery, remove red cable first. When connecting the battery, connect black cable first.

2. Extent of Delivery

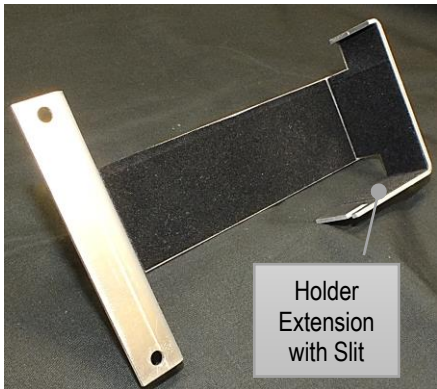
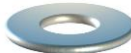



 <p>A: Battery Holder</p>	 <p>B: 2 x Washer Ø 10 mm</p>	 <p>C: 2 x Spacer Ø 8 mm</p>
	  <p>D: 2 x Adapter Cables (already connected to power supply cable of DIVING-PAM-II)</p>	
<p>E: 2 x Cable Ties</p>		

Figure 1: Extent of Delivery

3. Removing Battery Holder

See DIVING-PAM-II manual, Section “Battery Replacement” (https://www.walz.com/products/chl_p700/diving-pam-ii/downloads.html)

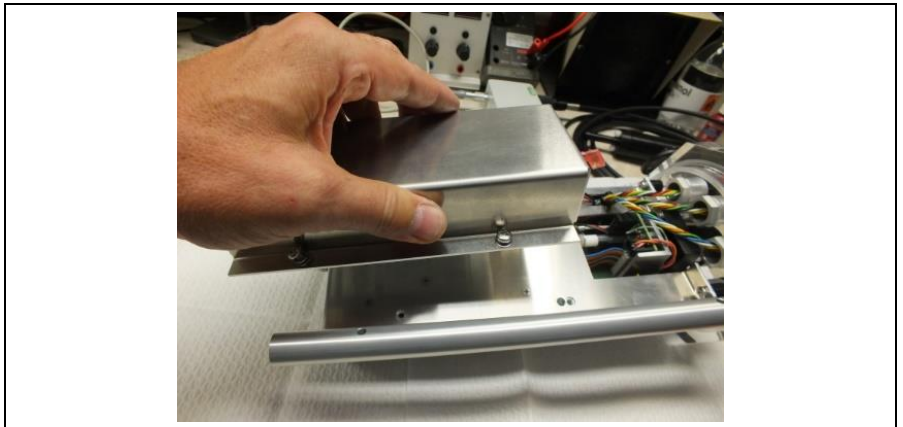


Figure 2: Remove Previous Holder

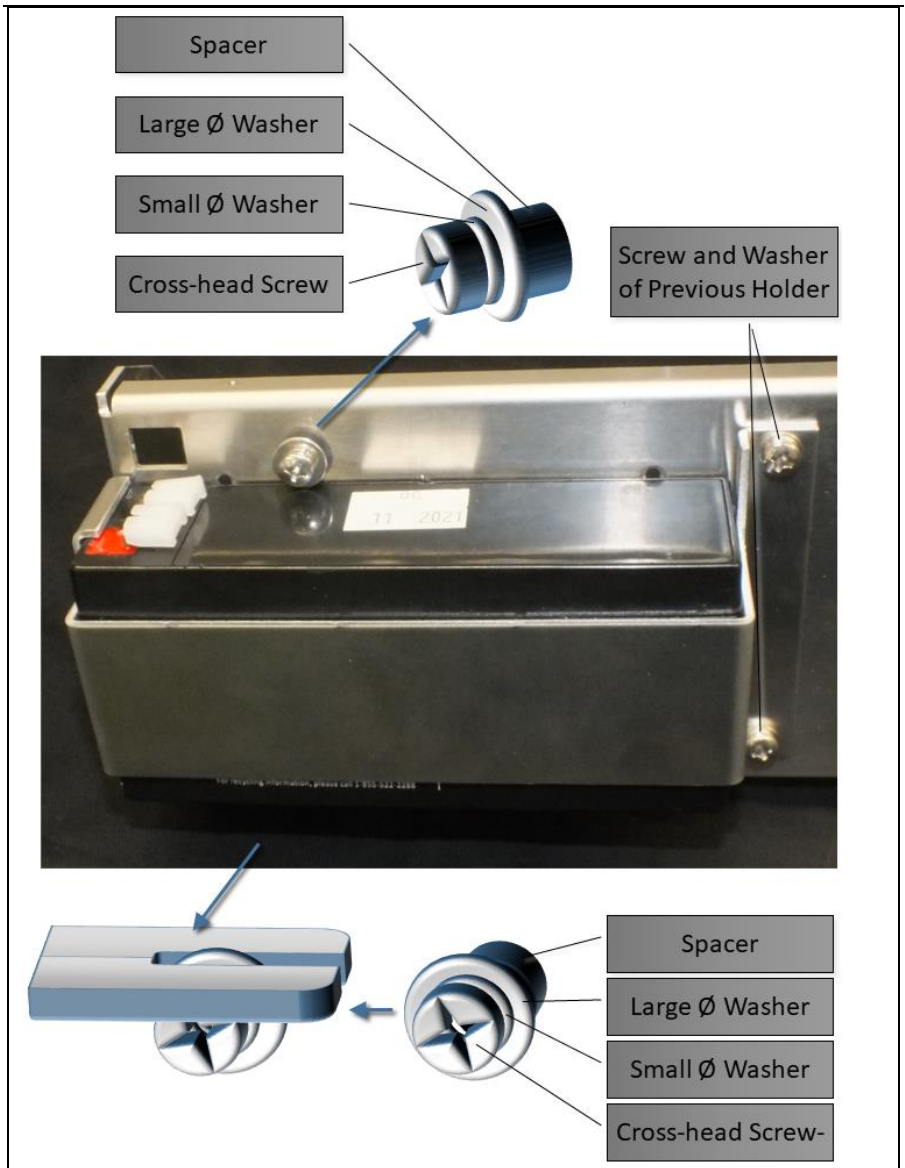


Figure 4: Mounting of Battery Holder

Not shown are the red and black charge cables connecting to the DIVING-PAM-II electronic board which are normally guided through the rectangular opening of the holder support.

5. Connecting the Battery

Firstly, connect adapter cables to power supply cables of DIVING-PAM-II (Figure 5). Then, remove protective caps from battery electrodes and connect adapter cables to battery electrodes. Finally, tie together cables and connections as shown in Figure 6.

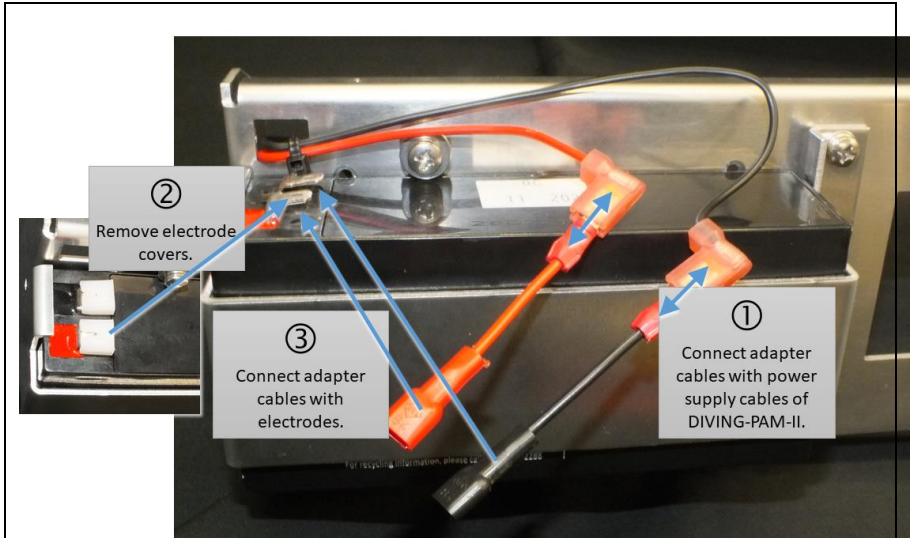


Figure 5: Connecting the Battery



Figure 6: Bundling using cable ties